
M3 LITE

BY **MOVO** 

USER MANUAL **GUIDE D'UTILISATION**

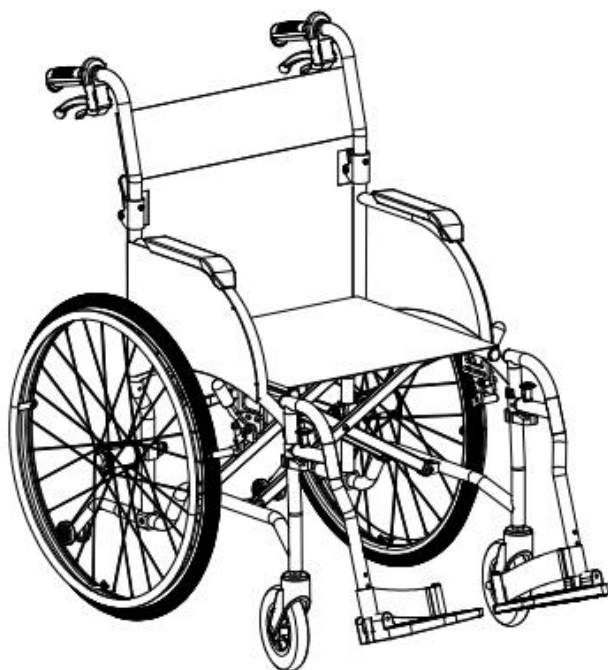MANUAL WHEELCHAIR **FAUTEUIL ROULANT**

Before operating this product, please read the user manual carefully and keep it for future reference.
Before operating this product, please read the user manual carefully and keep it for future reference.

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1. USER'S GUIDELINES

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Thank you for choosing our manual wheelchair. This product is made from specially processed magnesium and aluminum alloys, providing a strong yet lightweight design for greater comfort and ease of use.

Your safety is our priority. Please read this manual carefully before using the wheelchair, as it contains important instructions for safe operation.

This manual is an essential part of the product. Due to ongoing improvements, some illustrations may differ from your wheelchair. The manufacturer reserves the right to make changes without prior notice.

We hope this product enhances your daily life. We welcome your feedback on our wheelchair and after-sales service.

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2. SAFETY INFORMATION

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2.1 Unpacking Inspection

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Check that all components and accessories are included. If any parts are missing or damaged, contact your distributor or our company immediately.

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2.2 Symbol Descriptions

Description des symboles

Warning – Indicates a risk of serious injury or death if not used properly.

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Caution – Indicates a risk of minor injury or property damage if not used properly.

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Mandatory – Indicates required actions.

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Prohibited – Indicates actions that are not allowed.

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Maximum Load – Indicates the wheelchair's maximum weight capacity.

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Read Manual – Indicates that the manual must be read carefully.

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Pinch Hazard – Indicates a risk of finger injury during use.

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2.4 Intended Use

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This product is a mobility aid for individuals with reduced mobility, including persons with disabilities, patients, and older or physically weakened users.

It is designed for single-user use only, with a maximum weight of 120 kg.

Do not allow independent use by minors or by individuals with intellectual disabilities or who cannot operate the wheelchair safely.

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3.1 Product Features

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3.1.1 The wheelchair can be self-propelled or assisted by a caregiver or companion.

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3.1.2 This is a standard wheelchair. The backrest is not adjustable, the seat does not support position changes, and it is not intended for sports or bathing use.

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3.1.3 This product is intended for use in hospitals, care facilities, and home environments. Use on flat surfaces or gentle slopes. Recommended for dry surfaces free of water or ice.

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3.2 Main Structural Components

Composants principaux

The wheelchair consists of the frame, armrests, seat, wheels, leg rests, and other components.

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
3.2.1 Component Description (Figure 3-1)

3.2.1 Description des composants (figure 3-1)

- (1) Backrest upholstery
- (2) Push handles
- (3) Brake levers
- (4) Backrest folding mechanism
- (5) Armrests
- (6) Rear wheels
- (7) Anti-tip wheels
- (8) Wheel locks
- (9) Front wheels
- (10) Footplates
- (11) Leg rests
- (12) Seat upholstery

Figure 3-1

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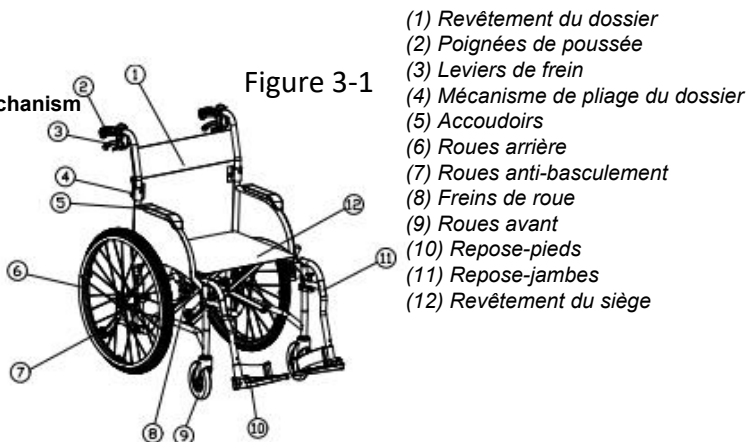INSTALLATION ET PLIAGE

mevo

4.1 Unfolding and Folding the Wheelchair

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4.1.1 Unfolding

4.1.1.1 Ensure the brakes are unlocked. Stand at the front or rear, grasp both armrests, and pull outward firmly. (Figure 4-1 (next page))

4.1.1.2 Press down on both sides of the seat until fully opened. (Figure 4-2 (next page))

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Figure 4-1

Figure 4-2

WARNING! Do not pull the armrests forcefully when unfolding, as this may damage the wheelchair. When pressing down on the seat, keep fingers clear of the seat tubes to avoid injury.

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4.1.2 Folding

Fold the footplates upward (C).

Lift the seat from the center using both sides of the seat upholstery (A).

Then push the armrests inward (B) to close the wheelchair. (Figure 4-3)

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4.2 Opening and Folding the Backrest

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4.1.1 Unfolding

Hold the push handles and lift them upward. Ensure that the locking pin on the backrest folding mechanism is fully engaged and visible. (Figures 4-4 and 4-5)

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Figure 4-5

WARNING! Avoid pinching fingers when folding or unfolding.

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4.2.2 Folding

Hold a push handle with one hand. With the other, press the backrest folding lever to release the locking pin. (Figure 4-6)

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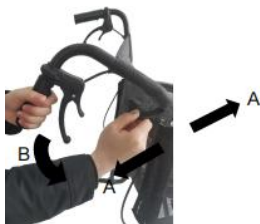
Figure 4-6

4.3 Installation and Removal of Leg Rests

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4.3.1 Installation

4.3.1.1 Insert the leg rest pivot (①) into the front frame tube (②) from the outside.


4.3.1.2 Rotate the leg rest inward until the latch pin (⑤) engages the mounting block (④). A “click” confirms proper installation.

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- ① Leg rest pivot
- ② Front frame tube
- ③ Release lever
- ④ Mounting block
- ⑤ Latch pin

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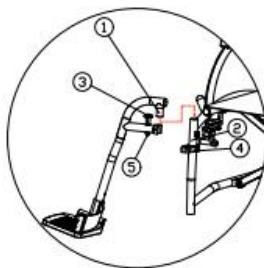4.3.2 Removal

4.3.2.1 Lift the release lever (③) to disengage the latch pin (⑤) from the mounting block (④).


4.3.2.2 Rotate the leg rest outward, then lift it off to remove.

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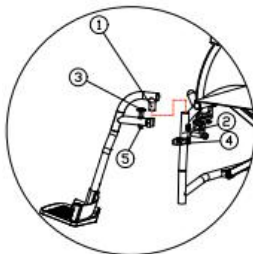4.3.3 Footplate Adjustment



4.3.3.1 Ensure at least 50 mm clearance between the footplates and the ground, and adjust both sides to the same height.

4.3.3.2 The footplates can be adjusted to five positions using the provided tool.

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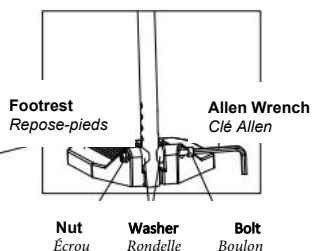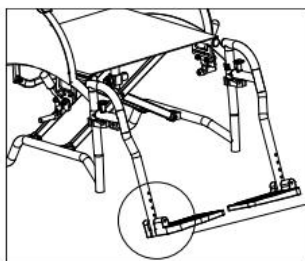4.4 Rear Wheel Installation and Removal

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
4.4.1 Installation

Press the quick-release button (①), insert the axle through the rear wheel (②) into the rear frame tube (③), and push until the locking balls fully extend.

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- ① Quick-release axle
- ② Rear wheel
- ③ Rear frame tube

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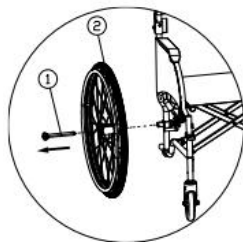
4.4.2 Removal

Press the button on the quick-release axle (①) and pull the rear wheel assembly (②) outward to remove it.

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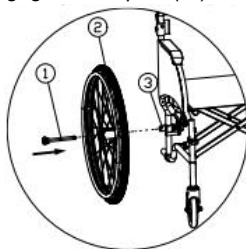4.5 Anti-Tip Wheel Installation and Removal

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
4.5.1 Installation

Press the release balls on the anti-tip wheel assembly (①), then insert it into the rear frame tube (②) until the balls fully extend.

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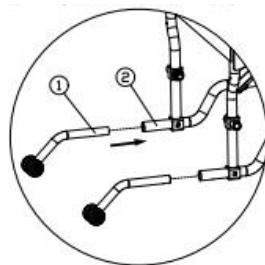- ① Anti-tip wheel assembly
- ② Rear frame tube

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
4.5.2 Removal

Press the release balls and pull the anti-tip wheel outward to remove.

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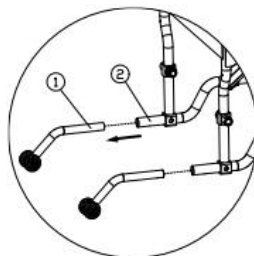4.6 Seat and Backrest Cushion Installation and Removal

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Attach the seat and backrest cushions using the hook-and-loop fasteners. Press firmly to secure. To remove, detach the fasteners.

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4.7 Transport and Handling

Transport et manutention

4.7.1 Move the wheelchair by pushing it, allowing the wheels to roll freely.

4.7.2 The wheelchair can be folded and lifted by the seat tubes.

4.7.3 For transport in a vehicle, fold the wheelchair and remove the rear wheels, anti-tip wheels, and footrests before placing it in the trunk.

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5.1 Self-Operation

Utilisation autonome

5.1.1 Propel the wheelchair by turning the handrims.

5.1.2 To go up a step, stop and position the front wheels at 90°. Lean forward and push firmly. Do not attempt steps over 30 mm alone.

5.1.3 On slopes, lean forward when going up. When going down, stay forward and control speed using the handrims. Do not attempt steep or long slopes alone.


5.1.4 Cross grooves or gaps at a right angle. Do not cross if too wide or deep.

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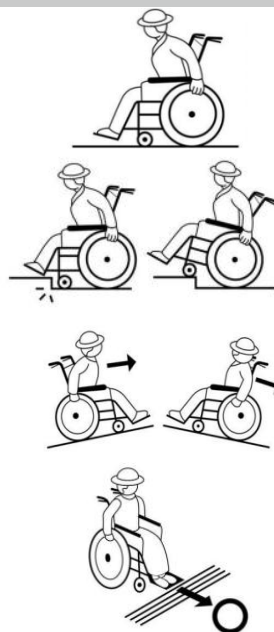5.2.1 On sidewalks, hold the handles firmly with both hands and push forward.

5.2.2 To go up a step, stop in front. Hold the handles, press down, and use your foot on the tipping bar to lift the front wheels while pushing forward. Adjust anti-tip wheels if equipped.

5.2.3 On slopes, lean slightly forward and push uphill. When descending, maintain this posture and move backward to control speed.

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
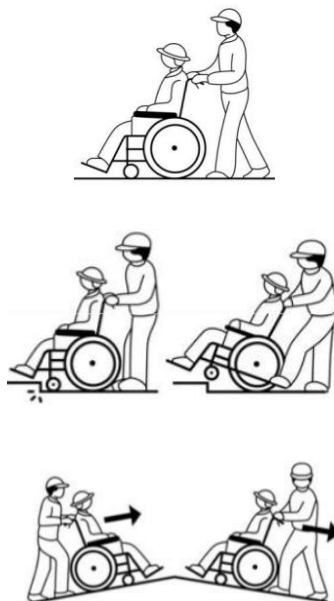
5.3 Manual Brake Operation

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While seated, pull both brake levers backward to engage. Push forward to release.

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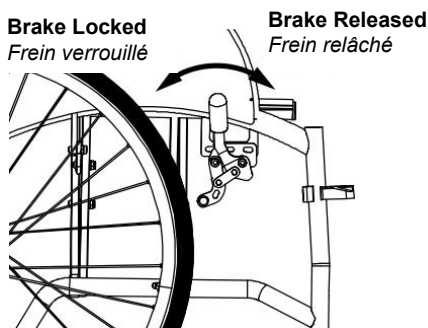5.4.1 Pull the brake lever up to engage. Release to disengage.

5.4.2 Lift lightly to slow the wheelchair.

5.4.3 Lift firmly to apply the parking (lock) brake.

5.4.4 To release the parking brake, lift the lock handle.

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5.5.1 Self-Operation

Utilisation autonome

- 5.5.1.1** Read the manual. Assemble and adjust the footrests, rear wheels, and anti-tip wheels.
- 5.5.1.2** Install the seat and backrest cushions.
- 5.5.1.3** Check that the brakes and backrest lock are secure.
- 5.5.1.4** Sit and adjust your position and balance.
- 5.5.1.5** With assistance, practice in a flat, open area.
- 5.5.1.6** Push the handrims forward to move ahead.
- 5.5.1.7** Pull the handrims backward to reverse.
- 5.5.1.8** To turn left, push the right handrim more; to turn right, push the left more.
- 5.5.1.9** Practice braking.
- 5.5.1.10** Practice transfers in and out of the wheelchair.

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5.5.2 Caregiver Operation

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- 5.5.2.1** Read the manual. Assemble and adjust the footrests, rear wheels, and anti-tip wheels.
- 5.5.2.2** Install the seat and backrest cushions.
- 5.5.2.3** Check that the brakes and backrest lock are secure.
- 5.5.2.4** Ensure the user is seated properly and balanced.
- 5.5.2.5** Practice pushing forward and backward and using the brakes to stop and hold the wheelchair.

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5.6 Transfers In and Out of the Wheelchair

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5.6.1 Position the wheelchair close to the transfer surface (e.g., bed or sofa).

5.6.2 Engage both brakes.

5.6.3 Flip up the footplates. Do not stand on them.

5.6.4 Transfer from the side.

5.6.5 Users with sufficient strength may stand or sit to transfer.

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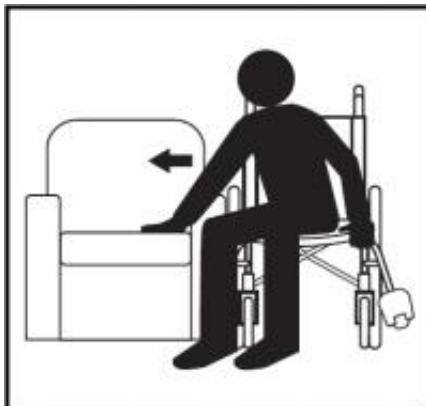
Note: Use assistance if needed.

Note: Always engage the brakes during transfers to prevent tipping.

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6.1 Cleaning

Nettoyage

6.1.1 Clean the seat and backrest cushions weekly with water and mild detergent.

6.1.2 Wipe other surfaces weekly with a clean cloth or disinfectant wipe.

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6.2 Disinfection

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6.2.1 Disinfect armrests, push handles, and similar parts weekly using a cloth or wipe with 75% medical alcohol. Dry immediately.

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Caution:

- a) Do not use high-pressure cleaners or spray water on the wheelchair.
- b) Do not use high-concentration alcohol, ether, strong acids, chlorinated solvents, gasoline, kerosene, or other harsh or abrasive cleaners.
- c) Ensure all parts are fully dry before use.

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- a) Check regularly for loose parts or bolts; tighten if needed.
 - b) Lubricate moving parts weekly.
 - c) Wipe off moisture and dirt after use to prevent rust.
 - d) Keep seat and backrest cushions clean.
 - e) Store in a dry, well-ventilated, non-corrosive place.
 - f) Ensure the frame folds smoothly.
 - g) Check that user and caregiver brakes work properly.
 - h) Ensure wheels rotate smoothly without wobble.
-
- a) *Vérifiez régulièrement si des pièces ou des boulons sont desserrés et resserrez au besoin.*
 - b) *Lubrifiez les pièces mobiles chaque semaine.*
 - c) *Essuyez l'humidité et la saleté après usage pour prévenir la rouille.*
 - d) *Gardez les coussins de siège et de dossier propres.*
 - e) *Entreposez dans un endroit sec, bien ventilé et non corrosif.*
 - f) *Assurez-vous que le cadre se plie facilement.*
 - g) *Vérifiez le bon fonctionnement des freins utilisateur et du proche aidant.*
 - h) *Assurez-vous que les roues tournent librement, sans jeu.*

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6.4 Reuse

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After long-term storage, clean, disinfect, and service the wheelchair per Sections 6.2 and 6.3 before reuse.

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6.5 Waste Disposal

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This symbol means the product must not be disposed of with household waste. At end of life, take it to an authorized collection or recycling facility. Proper disposal helps conserve resources and protect the environment. Components should be recycled per local regulations by authorized facilities. Plastic parts, the frame, and tubes can be recycled at end of life.

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7.1 Warranty and Replacement

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WARRANTY: 1 year from purchase. For products bought from authorized dealers or the company and used as instructed, defective parts due to materials or manufacturing will be repaired or replaced. Service must be performed at authorized centers.

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From manufacture date:

Frame: 3-year warranty

Consumables (e.g., tires, handgrips, seat and backrest cushions) are not covered; normal wear is the user's responsibility

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Warranty does not apply in the following cases:

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8. TRANSPORTATION AND STORAGE TRANSPORT ET ENTREPOSAGE

8.1 Transportation and Storage Conditions

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